


# DUCTILITY MACHINE

## ASTM D 113, D 5976 - IP 32 - ISO 1208

### EN 13398, 13589, 13703

The apparatus consists of a constant temperature bath capable to work in the range from -10 to 40 °C (thanks to a built-in cooling compressor) with an electrically powered tractor that pulls apart at a determined speed the two ends of two briquet specimens: the elongation is measured by means of a pointer on the tractor and a graduated scale fixed to the cover.

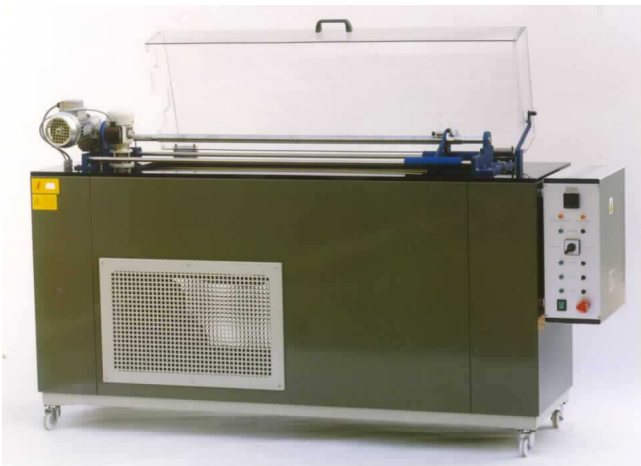

Apparatus

- Enamel finished steel case, floor mounted with swivelling castors.
- Stainless steel bath thermally insulated: complete with drain valve.
- Chromium plated brass pump: it provides an adequate stirring of the liquid in the bath.
- Stainless steel heater.
- Microprocessor thermoregulator with PID action: built-in digital display 0.1°C accuracy. Probe: Pt100 RTD. Working range: from -10 to 40°C.
- Safety device against overheating and low level with alarm lamp.
- Built-in hermetic cooling compressor: CFC/HCFC-free gas.
- Acrylic hinged cover: it helps to maintain a constant temperature in the bath providing a clear view of the inside.
- Traction carriage with two-speed motor for 10 and 50 mm/min: the carriage supports 2 briquet molds. Automatic stop safety switch that stops the carriage when the maximum elongation has been reached. Traction length: 50, 100 or 150 cm.
- White screen fixed to the bottom of the bath: it permits a more clear vision of the elongated specimen during the test.
- Easy access control box containing all the electronics: anodized aluminium control panel with english written indications.
- English written user manual with installation instructions.
- For 220 V/50 Hz connections: 1800 W power consumption.
- Dimensions (l x w x h): mm 1800 x 450 x 1300 approx (AD0113-200) mm 2400 x 450 x 1300 approx (AD0113-250). Weight: 250 kg approx.
- CE marked.

AD0113-100    Apparatus, 50 cm traction length  
 AD0113-110    Apparatus, 100 cm traction length  
 AD0113-120    Apparatus, 150 cm traction length

#### ACCESSORIES

AD0113-A00    Mold acc. to ASTM D 113 - IP 32 - IS 1208  
 AD0113-A01    Brass plate for mold assembling  
 AD5976-A00    Mold acc. to ASTM D 5976  
 CAL001        PT100 simulator  
 CAL003        Official Certificate for Pt100 simulator

Specifications may vary without notice.

**The apparatus includes the items listed aside the picture, accessories etc. should be purchased separately.**

#### Agence Nord:

ZA Object'Ifs Sud - Lot A3  
 6 Allée Emilie du Châtelet  
 14123 Ifs  
 tél : 02.31.34.50.74  
 fax : 02.31.34.55.17

#### Agence Sud:

Bât Le Venango. 392 Rue Jean Dausset  
 AGROPARC - BP11575  
 84916 Avignon Cédex 9  
 tél : 04.90.27.17.95 fax : 04.90.27.17.52



#### Agence Est:

Parc Club des Tanneries  
 2 Rue de la Faisanderie  
 67380 Lingolsheim  
 tél : 03.88.04.01.81  
 fax : 03.68.93.01.52

**www.deltalabo.fr**  
**info@deltalabo.fr**