


CONSISTOMETER

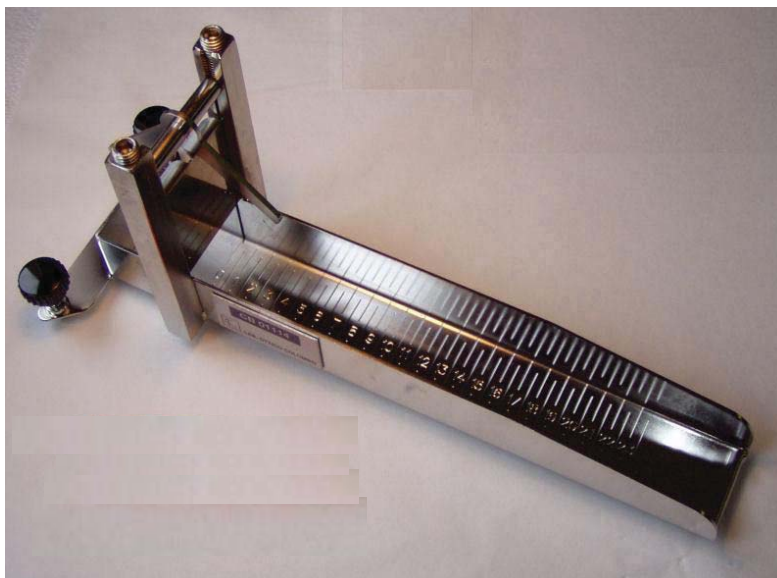

Bostwick Consistometer is useful to determine the consistency of cremogenates, concentrated juice and ketchup. Quick and certified method of international level. It gives an index by the rate of flow of a given volume of product, under constant conditions (sugar content and temperature) within a standard time.

TECHNICAL DATA:

- All Stainless-Steel Construction
- Reservoir Volume: 100ml
- Lane length: 23 cm
- Resolution: 0,5 cm
- Weight: 1200 g

Agence Nord:

Campus Effiscience
2 Rue Jean Perrin - Bât D
14460 Colombelles
tél : 02.31.34.50.74
fax : 02.31.34.55.17

Agence Sud:

Hôtel d'Entreprises de la Croix Rouge
Lot A4 - 10 Av de la Croix Rouge
84000 Avignon
tél : 04.90.27.17.95 fax : 04.90.27.17.52



Agence Est:

Parc Club des Tanneries
2 Rue de la Faisanderie
67380 Lingolsheim
tél /fax : 03.88.04.01.81

www.deltalabo.fr
info@deltalabo.fr