



Agence Nord:
ZA Object'ifs Sud - Lot A3
6 Allée Emilie du Châtelet
14123 Ifs
tél : 02.31.34.50.74
fax : 02.31.34.55.17



Agence Est:
Parc Club des Tanneries
2 Rue de la Faisanderie
67380 Lingolsheim
tél : 03.88.04.01.81
fax : 03.68.93.01.52

Agence Sud:
Bât Le Venango, 392 Rue Jean Dausset
AGROPARC - BP11575
84916 Avignon Cédex 9
tél : 04.90.27.17.95 fax : 04.90.27.17.52

www.deltalabo.fr
info@deltalabo.fr

Eco Titrator

Don't take risks in your
quality control

PEOPLE
YOU
CAN
TRUST

Metrohm
means ...
Spectroscopy!



 **Metrohm**

We give you peace of mind

When something goes wrong in quality control, the reputation of your brand is at stake. Better to exclude any risks in the first place.

With the new Metrohm Eco Titrator this is easy. The Eco Titrator makes titration so simple, safe, and reliable that everybody can perform titrations – and get the results right!

A SIMPLE AND RELIABLE ROUTINE

You can **save up to 120 titration methods** on your Eco Titrator. Simply select the required method, enter the sample size, and press the start button – that's it. The Eco Titrator performs the analysis fully automatically.

- Touch sensitive user interface
- Routine and expert user mode to protect your methods and prevent user errors
- Button to directly control the speed of the stirrer

SUPERIOR PRECISION

The **high resolution of the Eco Titrator's buret** allows precise dosing of volumes down to the micro liter range. This makes it easy to comply with the requirements of DIN EN ISO 8655-6.

- Flexible buret volume (5, 10, 20, or 50 mL)
- Straightforward, software guided installation and buret exchange



MULTILINGUAL

User-dialog in 13 languages: DE, EN, CN, FR, IT, JP, KR, PL, PT, RU, ES, TH, VN

COMPACT AND ROBUST

Great emphasis was put on the **compact and robust design** of the Eco Titrator. The complete device is easily moved from one place to another, so you always have it where you need it – be it in the laboratory or right next to the production line.

- Integrated magnetic stirrer
- Safe stand for the reagent bottle
- Straightforward installation and transportation

GLP COMPLIANT DOCUMENTATION

Using the USB port of the device, you can **connect a printer to your Eco Titrator** or use this port to store your methods and/or results on a USB memory stick.

- GLP compliant reports for immediate printing or storing on USB memory stick
- Archiving and easy method transfer using a USB memory stick.

InstantUp

With this function, you can put your Eco Titrator into a sleep mode (system off) when you are not using it. It is immediately ready for use again without you having to restart it completely.

RELIABLE RESULTS

The Eco Titrator features a **highly precise measuring** input. In combination with our electrodes «Made in Switzerland» this provides for measuring results of highest accuracy and reliability.

Eco Titrator – a true all-rounder

Whichever your industry is – the Eco Titrator is the right device for you. Combined with the right Metrohm sensor, you can perform all standard applications in potentiometric titration.

Connect a Custom Q3X thermal printer or a TX-900MH impact printer to the Eco Titrator's USB port and print out GLP compliant reports right in the laboratory.



Compact design:
just 28,6 cm (w) x 28,6 cm (d) x 22,0 cm (h)



WATER ANALYSIS

- Carbonate and calcium / magnesium hardness
- Chloride
- Sulfate
- Permanganate index




FOOD & BEVERAGES

- Acidity
- Chloride
- Vitamin C
- Iodine and peroxide number of fats and oil



PETRO-CHEMISTRY

- TAN / TBN
- Mercaptanes
- Chloride
- Bromine number

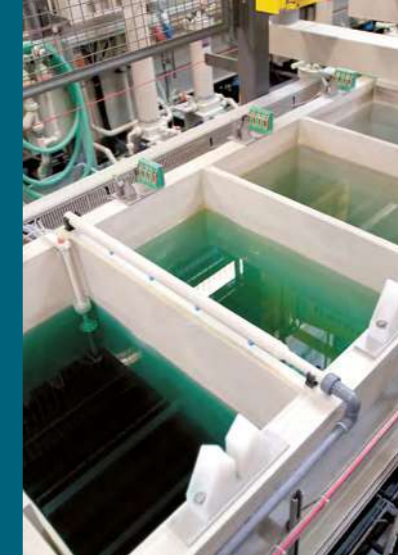

PLATING / SURFACE FINISHING

- TAN
- Metal content
- Chloride
- Analysis of degreasing / etching baths



PHOTOMETRIC TITRATIONS (with color indicator)

- Water hardness
- Metals in cement
- p- and m-value of water



SURFACTANT ANALYSIS

- Anionic / cationic surfactants
- Non ionic surfactants

Producing quality results can be so simple

Once your applications are saved on the device, working with the Eco Titrator could not be any easier: Press a few keys and what you get is an accurate, reproducible result for your sample. Depending on your requirements, save a GLP compliant report as a PDF file or send it to printer.



→ Select method



→ Enter the sample data and start the titration



→ Obtain a GLP compliant report



Local service and support – worldwide

With the Eco Titrator, you benefit from **quality «Made in Switzerland»** as well as professional service and support from our local sales organizations – wherever you are.

- Professional application support
- Technical service on site
- Certified Dosing Test



ORDERING INFORMATION

Instrument	
2.1008.0010	Eco Titrator
Accessories	
6.03002.150	Cylinder unit Eco 5 mL
6.03002.210	Cylinder unit Eco 10 mL
6.03002.220	Cylinder unit Eco 20 mL
6.03002.250	Cylinder unit Eco 50 mL
2.142.0100	Custom Q3X thermal printer
2.143.0200	Impact printer TX-900MH
Application specific packages	
These application specific packages include all required accessories for the dedicated applications. Unpack and go...	
2.1008.1010	Eco Titrator Acid / Base Complete package for aqueous pH titrations. Includes Eco Titrator, Cylinder unit 20 mL, combined pH electrode 6.0262.100 LL-Ecotrode plus and 6.2308.020 Electrolyte KCl 3 mol/L
2.1008.2010	Eco Titrator Salt Complete package for chloride titrations with Titrant AgNO ₃ . Includes Eco Titrator, Cylinder unit 10 mL, and 6.00430.100 Ag Titrode.
2.1008.3010	Eco Titrator Oil Complete package for non-aqueous pH titrations (e.g. acid number / base number of petrochemical products). Includes Eco Titrator, Cylinder unit 20 mL, and 6.0229.010 Solvotrode easyClean
2.1008.4010	Eco Titrator Redox Complete Package for Redox titrations (e.g. Iodometry, Cerimetry, Manganometry, etc.). Includes Eco Titrator, Cylinder unit 20 mL, and 6.0431.100 Pt Titrode
2.1008.5010	Eco Titrator Biogas Complete package for determining the FOS/TAC value for monitoring your biogas reactor. Includes Eco Titrator, Cylinder unit 20 mL, 6.00221.600 Ecotrode Gel with Pt1000, and all buffer solutions and electrolytes for calibration and storage of the electrodes.

Agence Nord:

ZA Object'Ifs Sud - Lot A3
6 Allée Emilie du Châtelet
14123 Ifs
tél : 02.31.34.50.74
fax : 02.31.34.55.17

Agence Sud:

Bât Le Venango. 392 Rue Jean Dausset
AGROPARC - BP11575
84916 Avignon Cédex 9
tél : 04.90.27.17.95 fax : 04.90.27.17.52



Agence Est:

Parc Club des Tanneries
2 Rue de la Faisanderie
67380 Lingolsheim
tél : 03.88.04.01.81
fax : 03.68.93.01.52

www.deltalabo.fr
info@deltalabo.fr